

### 3.13: Titrations

A **titration** is a volumetric technique in which a solution of one reactant (the titrant) is added to a solution of a second reactant (the "analyte") until the equivalence point is reached. The equivalence point is the point at which titrant has been added in exactly the right quantity to react stoichiometrically with the analyte (when moles of titrant = moles of analyte). If either the titrant or analyte is colored, the equivalence point is evident from the disappearance of color as the reactants are consumed. Otherwise, an **indicator** may be added which has an "endpoint" (changes color) at the equivalence point, or the equivalence point may be determined from a titration curve. The amount of added titrant is determined from its concentration and volume:

$$n \text{ (mol)} = C \text{ (mol/L)} * V \text{ (L)}$$

and the amount of titrant can be used in the usual stoichiometric calculation to determine the amount of analyte.

The titration process can be observed in the video below.



A measured volume of the solution to be titrated, in this case, colorless aqueous acetic acid,  $\text{CH}_3\text{COOH}(aq)$  is placed in a beaker. The colorless sodium hydroxide  $\text{NaOH}(aq)$ , which is the **titrant**, is added carefully by means of a buret. The volume of titrant added can then be determined by reading the level of liquid in the buret before and after titration. This reading can usually be estimated to the nearest hundredth of a milliliter, so precise additions of titrant can be made rapidly.



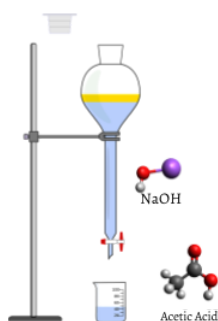


Figure 3.13.1: The titration setup initially, before titrant ( $\text{NaOH}$ ) has been added.  $\text{NaOH}$  is held in the buret, which is positioned above the beaker of acetic acid. Titrant ( $\text{NaOH}$ ) is added until it neutralizes all of the analyte (acetic acid). This is called the equivalence point. Note: Unlike the picture, both substances are actually clear but are blue for visibility purposes in the picture.

As the first few milliliters of titrant flow into the flask, some indicator briefly changes to pink, but returns to colorless rapidly. This is due to a large excess of acetic acid. The limiting reagent  $\text{NaOH}$  is entirely consumed.

The added indicator changes to pink when the titration is complete, indicating that all of the aqueous acetic acid has been consumed by  $\text{NaOH}(aq)$ . The reaction which occurs is

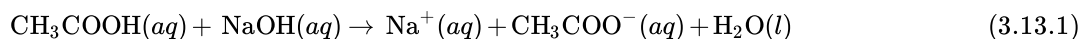

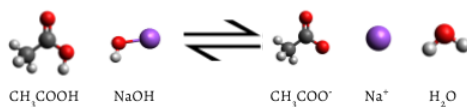

Reaction of acetic acid and sodium hydroxide to give acetate ion, sodium ion and water. The reaction is shown in terms of stick and ball diagram of each species.

Eventually, all the acetic acid is consumed. Addition of even a fraction of a drop of titrant produces a lasting pink color due to unreacted NaOH in the flask. The color change that occurs at the **endpoint** of the indicator signals that all the acetic acid has been consumed, so we have reached the **equivalence point** of the titration. If slightly more NaOH solution were added, there would be an excess and the color of the solution in the flask would get much darker. The endpoint appears suddenly, and care must be taken not to overshoot the endpoint.

After the titration has reached the endpoint, a final volume is read from the buret. Using the initial and final reading, the volume added can be determined quite precisely:

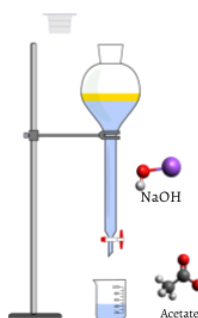


Figure 3.13.2 The figure above shows a completed titration, where the equivalence point has been reached. NaOH (the titrant) has neutralized all of the Acetic Acid, leaving Acetate in the beaker. At this point, the moles of NaOH added is equivalent to the moles of acetic acid initially in the beaker.

The object of a titration is always to add just the amount of titrant needed to consume *exactly* the amount of substance being titrated. In the NaOH—CH<sub>3</sub>COOH reaction Eq. 3.13.1, the equivalence point occurs when an equal molar amount of NaOH has been added from the graduated cylinder for every mole of CH<sub>3</sub>COOH originally in the titration flask. That is, at the equivalence point the ratio of the amount of NaOH, added to the amount of CH<sub>3</sub>COOH consumed must equal the stoichiometric ratio

$$\frac{n_{\text{NaOH}}(\text{added from graduated cylinder})}{n_{\text{CH}_3\text{COOH}}(\text{initially in flask})} = S\left(\frac{\text{NaOH}}{\text{CH}_3\text{COOH}}\right)$$

$$= \frac{1 \text{ mol NaOH}}{1 \text{ mol CH}_3\text{COOH}}$$

#### ✓ Example 3.13.1 : Endpoint of Titration

What volume of 0.05386 M KMnO<sub>4</sub> would be needed to reach the endpoint when titrating 25.00 ml of 0.1272 M H<sub>2</sub>O<sub>2</sub>, given S(KMnO<sub>4</sub>/H<sub>2</sub>O<sub>2</sub>) = 2/5

#### Solution

At the equivalence point, the stoichiometric ratio will apply, and we can use it to calculate the amount of KMnO<sub>4</sub> which must be added:

$$n_{\text{KMnO}_4}(\text{added}) = n_{\text{H}_2\text{O}_2}(\text{in flask}) \times S\left(\frac{\text{KMnO}_4}{\text{H}_2\text{O}_2}\right)$$

The amount of H<sub>2</sub>O<sub>2</sub> is obtained from the volume and concentration:

$$n_{\text{H}_2\text{O}_2}(\text{in flask}) = 25.00 \text{ cm}^3 \times 0.1272 \frac{\text{mmol}}{\text{cm}^3} = 3.180 \text{ mmol H}_2\text{O}_2$$

Then

$$\begin{aligned} n_{\text{KMnO}_4}(\text{added}) &= 3.180 \text{ mmol H}_2\text{O}_2 \times \frac{2 \text{ mol KMnO}_4}{5 \text{ mol H}_2\text{O}_2} \times \frac{10^{-3}}{10^{-3}} \\ &= 3.180 \text{ mmol H}_2\text{O}_2 \times \frac{2 \text{ mmol KMnO}_4}{5 \text{ mmol H}_2\text{O}_2} \end{aligned}$$

$$= 1.272 \text{ mmol KMnO}_4$$

To obtain  $V_{\text{KMnO}_4(aq)}$  we use the concentration as a conversion factor:

$$V_{\text{KMnO}_4(aq)} = 1.272 \text{ mmol KMnO}_4 \times \frac{1 \text{ cm}^3}{5.386 \times 10^{-2} \text{ mmol KMnO}_4}$$


$$= 23.62 \text{ cm}^3$$

Note that overtitrating [adding more than  $23.62 \text{ cm}^3$  of  $\text{KMnO}_4(aq)$ ] would involve an excess (more than 1.272 mmol) of  $\text{KMnO}_4$ .

Titration is often used to determine the concentration of a solution. In many cases it is not a simple matter to obtain a pure substance, weigh it accurately, and dissolve it in a volumetric flask as was done in [Example 1 of Solution Concentrations](#). NaOH, for example, combines rapidly with  $\text{H}_2\text{O}$  and  $\text{CO}_2$  from the air, and so even a freshly prepared sample of solid NaOH will not be pure. Its weight would change continuously as  $\text{CO}_2(g)$  and  $\text{H}_2\text{O}(g)$  were absorbed. Hydrogen chloride (HCl) is a gas at ordinary temperatures and pressures, making it very difficult to handle or weigh. Aqueous solutions of both of these substances must be **standardized**; that is, their concentrations must be determined by titration.

### ✓ Example 3.13.2: Concentration of Titrant

A sample of pure potassium hydrogen phthalate ( $\text{KHC}_8\text{H}_4\text{O}_4$ ) weighing 0.3421 g is dissolved in distilled water. Titration of the sample requires 27.03 ml  $\text{NaOH}(aq)$ . The titration reaction is

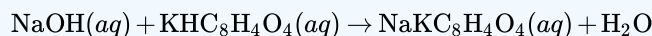

What is the concentration of  $\text{NaOH}(aq)$  ?

### Solution

To calculate concentration, we need to know the amount of NaOH and the volume of solution in which it is dissolved. The former quantity could be obtained via a stoichiometric ratio from the amount of  $\text{KHC}_8\text{H}_4\text{O}_4$ , and that amount can be obtained from the mass

$$\begin{aligned} m_{\text{KHC}_8\text{H}_4\text{O}_4} &\xrightarrow{M_{\text{KHC}_8\text{H}_4\text{O}_4}} n_{\text{KHC}_8\text{H}_4\text{O}_4} \xrightarrow{S(\text{NaOH}/\text{KHC}_8\text{H}_4\text{O}_4)} n_{\text{NaOH}} \\ n_{\text{NaOH}} &= 0.3421 \text{ g} \times \frac{1 \text{ mol KHC}_8\text{H}_4\text{O}_4}{204.22 \text{ g}} \times \frac{1 \text{ mol NaOH}}{1 \text{ mol KHC}_8\text{H}_4\text{O}_4} \\ &= 1.674 \times 10^{-3} \text{ mol NaOH} = 1.675 \text{ mmol NaOH} \end{aligned}$$

The concentration is


$$c_{\text{NaOH}} = \frac{n_{\text{NaOH}}}{V} = \frac{1.675 \text{ mmol NaOH}}{27.03 \text{ cm}^3} = 0.06197 \text{ mmol cm}^{-3}$$

or 0.06197 M.

By far the most common use of titrations is in determining unknowns, that is, in determining the concentration or amount of substance in a sample about which we initially knew nothing. The next example involves an unknown that many persons encounter every day.

### ✓ Example 3.13.3 : Vitamin C

Vitamin C tablets contain ascorbic acid ( $\text{C}_6\text{H}_8\text{O}_6$ ) and a starch “filler” which holds them together. To determine how much vitamin C is present, a tablet can be dissolved in water and with sodium hydroxide solution,  $\text{NaOH}(aq)$ . The equation is

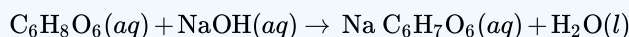

If titration of a dissolved vitamin C tablet requires  $16.85 \text{ cm}^3$  of  $0.1038 \text{ M NaOH}$ , how accurate is the claim on the label of the bottle that each tablet contains  $300 \text{ mg}$  of vitamin C?

### Solution

The known volume and concentration allow us to calculate the amount of  $\text{NaOH}(aq)$  which reacted with all the vitamin C. Using the stoichiometric ratio

$$S\left(\frac{\text{C}_6\text{H}_8\text{O}_6}{\text{NaOH}}\right) = \frac{1 \text{ mmol C}_6\text{H}_8\text{O}_6}{1 \text{ mmol NaOH}}$$

we can obtain the amount of  $\text{C}_6\text{H}_8\text{O}_6$ . The molar mass converts that amount to a mass which can be compared with the label. Schematically

$$V_{\text{NaOH}} \rightarrow c_{\text{NaOH}} n_{\text{NaOH}} \rightarrow S(\text{C}_6\text{H}_8\text{O}_6/\text{NaOH}) n_{\text{C}_6\text{H}_8\text{O}_6} \rightarrow M_{\text{C}_6\text{H}_8\text{O}_6} m_{\text{C}_6\text{H}_8\text{O}_6} \quad (3.13.2)$$

$$m_{\text{C}_6\text{H}_8\text{O}_6} = 16.85 \text{ cm}^3 \times \frac{0.1038 \text{ mmol NaOH}}{1 \text{ cm}^3} \times \frac{1 \text{ mmol C}_6\text{H}_8\text{O}_6}{1 \text{ mmol NaOH}} \times \frac{176.1 \text{ mg}}{\text{mmol C}_6\text{H}_8\text{O}_6} \quad (3.13.3)$$

$$= 308.0 \text{ mg} \quad (3.13.4)$$

Note that the molar mass of  $\text{C}_6\text{H}_8\text{O}_6$

$$\begin{aligned} \frac{176.1 \text{ g}}{1 \text{ mol C}_6\text{H}_8\text{O}_6} &= \frac{176.1 \text{ g}}{1 \text{ mol C}_6\text{H}_8\text{O}_6} \times \frac{10^{-3}}{10^{-3}} \\ &= \frac{176.1 \text{ g} \times 10^{-3}}{10^{-3} \text{ mol C}_6\text{H}_8\text{O}_6} = \frac{176.1 \text{ mg}}{1 \text{ mmol C}_6\text{H}_8\text{O}_6} \end{aligned}$$

can be expressed in milligrams per millimole as well as in grams per mole.

The  $308.0 \text{ mg}$  obtained in this example is in reasonably close agreement with the manufacturer’s claim of  $300 \text{ mg}$ . The tablets are stamped out by machines, not weighed individually, and so some variation is expected.

This page titled [3.13: Titrations](#) is shared under a [CC BY-NC-SA 4.0](#) license and was authored, remixed, and/or curated by [Ed Vitz, John W. Moore, Justin Shorb, Xavier Prat-Resina, Tim Wendorff, & Adam Hahn](#).